

# Model WIC 3 | Walk-in-chambers for stability testing according to ICH (Q1A)

The walk-in-chambers offer a wide range of features and benefits to perfectly meet requirements in stability testing. Perfect test conditions are maintained in a wide range of temperature from 10°C to 50°C and relative humidity of 20% to 90% RH.

## BENEFITS

- Stable and uniform environment
- Fully complies with Pharma regulations (ICH)Equipment
- Longevity - corrosion resistant
- Ease of operation
- Environmental friendly refrigerant with GWP value 547
- Maintenance possible with the least possible disruption to the climate conditioning
- Perfectly meets requirements in Pharma industry
- 3 years warranty



Model 3



Model 3

## IMPORTANT FEATURES

- Relative humidity accuracy of  $\pm 2,5\%$
- Temperature accuracy of  $\pm 1,5^\circ\text{C}$
- Duct work for proper distribution of air
- Highly insulated room panels with U value of only  $0,21 \text{ W/m}^2\text{K}$
- LED Lighting with motion detector
- High quality room elements, made in Germany
- Interior completely stainless steel
- Exterior service accessibility
- Room Elements: Galvanized sheet steel, polyester coated, white, similar to RAL 9010, interior surfaces in stainless steel, 4 inch / 100 mm PU foam U value:  $0,21\text{W/m}^2\text{K}$
- Floor: Floor overlay is chrome-nickel steel, pattern-rolled 5WL
- Surface load:  $50.000 \text{ N/m}^2$

## ORDERING INFORMATION


Interior volume [L]	Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
<b>Model WIC 3</b>					
24,516	230 V 1~ 50 Hz -20,0 A	fixed supply	Standard	WIC3-230V	9800-0013
			Standard	WIC3-208V	9800-0007
	208...230 V 1~ 60 Hz -20,0 A	fixed supply	with conditioner on ceiling	WIC3C-208V	9800-0008

## TECHNICAL DATA


Data			
Designation	WIC3-230V	WIC3-208V	WIC3C-208V
Article Number	9800-0013	9800-0007	9800-0008
Option model	Standard	Standard	with conditioner on ceiling
Performance Data Temperature			
Temperature range	{3299}...{3298}	{3299}...{3298}	{3299}...{3298}
Performance Data Climate			
Temperature range with humidity	10...50 °C	10...50 °C	10...50 °C
Humidity range	20...90 % RH	20...90 % RH	20...90 % RH
Temperature uniformity at 25°C and 60% RH	1,5 ± K	1,5 ± K	1,5 ± K
Temperature uniformity at 40°C and 75% RH	1,5 ± K	1,5 ± K	1,5 ± K
Temperature fluctuation at 25°C and 60% RH	0,5 ± K	0,5 ± K	0,5 ± K
Temperature fluctuation at 40°C and 75% RH	0,5 ± K	0,5 ± K	0,5 ± K
Electrical data			
Rated Voltage	230 V	208...230 V	208...230 V
Power frequency	50 Hz	60 Hz	60 Hz
Nominal power	2,7 kW	2,7 kW	2,7 kW
Unit fuse	20,0 A	20,0 A	20,0 A
Phase (Nominal voltage)	1~	1~	1~
Dimensions and weights			
Interior volume	24,516 L	24,516 L	24,516 L
Net weight of the unit (empty)	1,186 kg	1,186 kg	1,186 kg
Housing dimensions not incl. fittings and connections			
Width net	2,600 mm	2,600 mm	2,600 mm
Height net	2,470 mm	2,470 mm	2,470 mm
Depth net	4,700 mm	4,700 mm	4,700 mm
Internal Dimensions			
Interior width	2,400 mm	2,400 mm	2,400 mm
Interior height	2,270 mm	2,270 mm	2,270 mm

Interior depth	4,500 mm	4,500 mm	4,500 mm
Unit doors	1	1	1
<b>Environment-specific data</b>			
Energy consumption at 40°C and 75% RH	1,150 Wh/h	1,150 Wh/h	1,150 Wh/h

**DIMENSIONS Incl. fittings and connections [mm]**

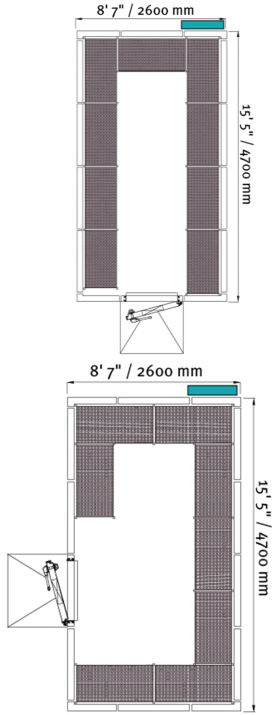

**DIAGRAMS**

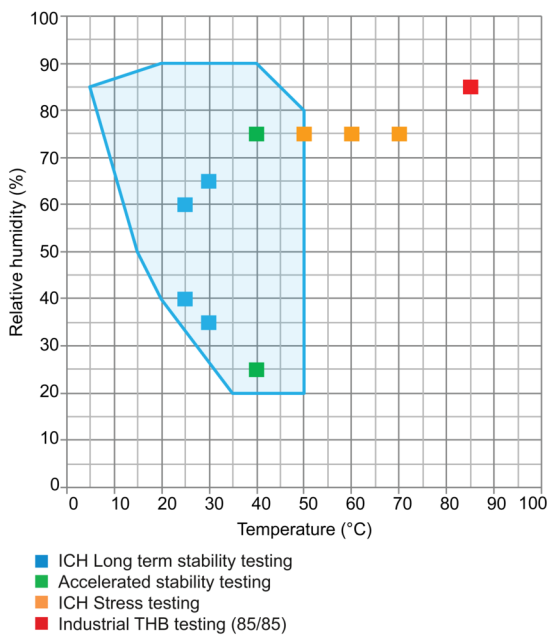

Climate chart

## ACCESSORIES

Designation	Description	WIC 3	*	Art.-No.
Measuring station	with 1x PT100 for humidity	•	-	8012-2579
Options/Accessories		•	-	8012-2461
Qualification documents	<p>IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature and humidity values</p> <p>Digital in PDF format</p> <p>Hard copy inside folder</p>	•	-	7057-0002
		•	-	7007-0002
Ramp for all WIC models	suitable for all WIC models	•	-	8012-2155
		•	-	8012-2239
Shelf set large	suitable for WIC 3, WIC 3 C	•	-	8012-2158
		•	-	8012-2242

## SERVICES

Designation	Description	*	Art.-No.
Access port with silicone plug	1 access port with silicone plug, 30 mm - position of your choice	-	DL10-1000
	1 access port with silicone plug, 50 mm - your choice of position	-	DL10-1050
	2 access ports with silicone plugs, 30 mm - your choice of position	-	DL10-1010
	2 access ports with silicone plugs, 50 mm - your choice of position	-	DL10-1060
	3 access ports with silicone plugs, 30 mm - your choice of position	-	DL10-1020
	3 access ports with silicone plugs, 50 mm - your choice of position	-	DL10-1070
	4 access ports with silicone plugs, 30 mm - your choice of position	-	DL10-1030
	4 access ports with silicone plugs, 50 mm - your choice of position	-	DL10-1080
	5 access ports with silicone plugs, 30 mm - your choice of position	-	DL10-1040
	5 access ports with silicone plugs, 50 mm - your choice of position	-	DL10-1090
<b>Calibration services</b>			
Temperature and humidity measurement, 13-1 measuring points	Temperature measurement with 13 temperature measuring points and 1 humidity measuring point in center of chamber with a pair of values specified by the user, including certificate	-	DL30-0113
<b>Installation services</b>			
Installation and Commissioning	Installation and Commissioning	-	DL10-0800
<b>Maintenance services</b>			

Designation	Description	*	Art.-No.
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	–	DL20-1000
<b>Validation services</b>			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder	–	DL40-0200

**BINDER GmbH**

Tuttlingen, Germany  
TEL +49 7462 2005 0  
info@binder-world.com  
www.binder-world.com

**BINDER Inc.**

Bohemia, NY, USA  
TEL +1 631 224 4340  
usa@binder-world.com  
www.binder-world.us

**BINDER Environmental Testing  
Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China  
TEL +86 21 685 808 25  
china@binder-world.com  
www.binder-world.com

**BINDER Asia Pacific (Hong Kong) Ltd.**

Kowloon, Hong Kong, P.R. China  
TEL +852 39070500  
asia@binder-world.com  
www.binder-world.com