

Model FDL 115 | Safety drying chambers for limited quantities of solvent

The BINDER FDL 115 with silicone and dust-free inner chamber, as well as symmetric airflow dries specimens containing solvents and thus this safety drying oven meets all EN 1539 requirements.

BENEFITS

- Safety concept meets the DIN EN 1539-2010-o8 standard
- Short drying time with high temperature accuracy and defined ventilation
- Silicone free and dust-free inner chamber
- Replaceable fresh air cartridges are easy to change and ensure constant fresh air quality



Model 115



Model 115

IMPORTANT FEATURES

- Temperature range: +10°C above ambient temperature to +300 °C
- APT.line™ preheating chamber technology
- Controller with time-segment and real-time programming
- 2 chrome-plated racks
- Fresh-air monitoring with audible and visual alarm and automatic deactivation of heating
- Replaceable filter cartridge, Class M6
- Class 2 independent adjustable temperature safety device (DIN 1288o) with visual alarm
- Computer interface: RS 422
- Standard door gasket made from solvent-resistant FKM for use at up to 200 °C
- Optional door gasket made from silicone with limited solvent resistance for use at up to 300 °C

ORDERING INFORMATION

Interior volume		Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
[L]	[L]					
Model FDL 115						
115	115	230 V 1~ 50/60 Hz -	CEE 7/7	Standard	FDL115-230V	9010-0292

TECHNICAL DATA

Data	
Designation	FDL115-230V
Article Number	9010-0292
Option model	Standard
Performance Data Temperature	
Temperature range	+10 °C above ambient temperature to 300 °C
Temperature uniformity at 150°C	2,5 ± K
Temperature fluctuation at 150°C	0,8 ± K
Heating up time to 150°C	20 min
Recovery time after door was opened for 30 s at 150°C	12 min
Air change data	
Air circulation (approx.)	20 x/min
Volumetric flow rate of exhaust air acc. to EN 1539 at 50 °C	400 ca. L/min
Solvent quantity	
Highest permitted solvent quantity (at T=180°C, M=100g/mol, U=40g/m ³ , K=0,5)	4,67 g
Electrical data	
Rated Voltage	230 V
Power frequency	50/60 Hz
Nominal power	2,9 kW
Phase (Nominal voltage)	1~
Measures	
Interior volume	115 L
Net weight of the unit (empty)	90 kg
Load per rack	30 kg
Permitted load	60 kg
Wall clearance back	160 mm
Wall clearance sidewise	100 mm
Housing dimensions not incl. fittings and connections	
Width net	830 mm
Height net	805 mm

Depth net	685 mm
-----------	--------

Internal Dimensions

Interior width	600 mm
----------------	--------

Interior height	435 mm
-----------------	--------

Interior depth	435 mm
----------------	--------

Environment-specific data

Sound-pressure level	57 dB(A)
----------------------	----------


Energy consumption at 15°C	1,200 Wh/h
----------------------------	------------

Fixtures

Number of shelves (std./max.)	2/5
-------------------------------	-----

All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10 %. The temperature data is determined in accordance to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. Technical data refers to 100 % fan speed. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS Incl. fittings and connections [mm]

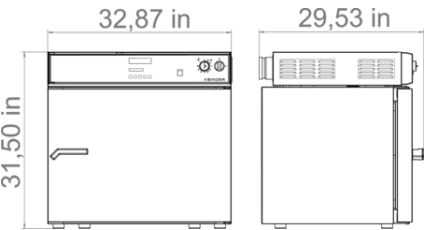

OPTIONS

Designation	Description	*	Art.-No.
Access port with silicone plug	back		
	10 mm	01	8012-1284
	left		
	10 mm	01	8012-1277
	30 mm	01	8012-1360
	right		
	10 mm	01	8012-1270
	30 mm	01	8012-1355
Calibration certificate, expanded	for temperature; for extending the measurement in center of chamber to include another test temperature	-	8012-1117

Designation	Description	*	Art.-No.
Calibration certificate, temperature	for temperature, measurement in center of chamber at specified temperature	-	8012-1136
	temperature measurement incl. certificate, 9 measuring points at specified temperature	-	8012-1553
	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	-	8012-1574
	temperature measurement incl. certificate and 27 measuring points at specified temperature	-	8012-1595
Coil coating door flap	door flap for coil coating tests	-	8012-0179
Door gasket	made from silicone, resistant to high temperatures for operating temperatures above 200°C. Note: The gasket wears out more quickly in temperatures above 250°C	-	8012-0499
Door lock	lockable door handle	-	8012-1248
Measurement of air exchange rate	in accordance with ASTM D5374, definition and protocol according to ambient temperature	-	8012-1196

ACCESSORIES

Designation	Description	*	Art.-No.
APT-COM™ 4	Multi Management Software APT-COM™		
	version 4, BASIC edition	–	9053-0039
	version 4, GLP edition	–	9053-0042
	version 4, PROFESSIONAL edition	–	9053-0040
Qualification documents	IQ/OQ/PQ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO ₂ , O ₂ – or pressure, depending on unit		
	Digital in PDF format	-	7057-0005
	Hard copy inside folder	-	7007-0005
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO ₂ , O ₂ , pressure, depending on unit		
	Digital in PDF format	-	7057-0001
	Hard copy inside folder	-	7007-0001
Rack	chrome plated	-	8012-2034
	stainless steel	-	8012-2048
Shelf, perforated	Stainless steel	-	8012-2161
Table on castors	stable cart, casters with locking brakes, dimensions: W 1,000 x D 800 x H 780 mm	-	9051-0018
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	8012-2250

SERVICES

Designation	Description	*	Art.-No.
Calibration services			
Measurement of air exchange rate	including certificate (in accordance with ASTM D5374)	–	DL33-0000
Temperature calibration	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate Travel to be quoted separate.	–	DL30-0101
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate Travel to be quoted separate.	–	DL30-0102
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	–	DL30-0109
Temperature measurement, 18 measuring points	Temperature measurement with 18 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 2 shelves are required.	–	DL30-0118
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	–	DL30-0127
Installation services			
Unit installation	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work) Travel to be quoted as separate.	–	DL10-0100
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming. Travel to be quoted as separate Limited to 3 hours.	–	DL10-0500
Maintenance services			
BRONZE 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate. 10% Discount on Spare Parts and Labor for any additional Services needed Wear + tear parts to be purchased or billed to customer. Travel to be quoted as separate.	–	DL20-0710
GOLD 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate of Temperature. 25% Discount on Spare Parts and Labor for any additional Services needed. One additional courtesy visit per year can be used for service visit or operator. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	–	DL20-0910
Preventive maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Checking a test temperature specified by the user in center of usable space without certificate. Wear + tear parts need to be purchased or billed to customer. Travel to be quoted as separate.	–	DL20-0200
SILVER 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate 15% Discount on Spare Parts and Labor for any additional Services needed One additional courtesy visit per year can be used for service visit or operator training. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	–	DL20-0810
Validation services			
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder. Travel to be quoted separate.	–	DL40-0100
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder. Travel to be quoted separate.	–	DL44-0500

Designation	Description	*	Art.-No.
Warranty service			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	–	DL50-0010

BINDER GmbH

Tuttlingen, Germany
TEL +49 7462 2005 0
FAX +49 7462 2005 100
info@binder-world.com
www.binder-world.com

BINDER Inc.

Bohemia, NY, USA
TEL +1 631 224 4340
FAX +1 631 224 4354
usa@binder-world.com
www.binder-world.us

**BINDER Environmental Testing
Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China
TEL +86 21 685 808 25
FAX +86 21 685 808 29
china@binder-world.com
www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.

Kowloon, Hong Kong, P.R. China
TEL +852 39070500
FAX +852 39070507
asia@binder-world.com
www.binder-world.com