

## Model BD-S 56 | Standard-Incubators with gravity convection

The new products in the Solid.Line impress with their tried-and-tested BINDER quality and reliability. They are suited to incubation. The attributes of the Solid.Line are geared towards the application in question and therefore used in research and quality assurance.

### BENEFITS

- Easy handling
- Accurate temperature control
- Simple and ergonomic door opening



Model 56

### IMPORTANT FEATURES

- Temperature range: ambient temperature +5 °C to 70 °C
- APT.line™ preheating chamber technology
- Gravity convection
- Adjustable exhaust air flap
- Controller with timer function
- Inner door made of tempered safety glass
- 1 chrome-plated rack, incl. shelf supports
- Class 3.1 integrated independent temperature safety device (DIN 12880) with visual alarm

### ORDERING INFORMATION


Interior volume	Power supply - unit fuse	Plug*	Version	Model version	Art.-No.
[L]	[L]				
<b>Model BD-S 56</b>					
62	62	230 V 1~ 50/60 Hz -6,3 A	CEE 7/7	Standard	BDS056-230V
		120 V 1~ 60 Hz -12,5 A	NEMA 5-15	Standard	BDS056UL-120V

## TECHNICAL DATA

Data		
Designation	BDS056-230V	BDS056UL-120V
Article Number	9090-0016	9090-0017
Option model	Standard	Standard
Performance Data Temperature		
Temperature range	+5 °C above ambient temperature to 70 °C	+5 °C above ambient temperature to 70 °C
Temperature uniformity at 37°C	0.5 ± K	0.5 ± K
Temperature fluctuation at 37°C	0.3 ± K	0.3 ± K
Heating up time to 37°C	230 min	230 min
Recovery time after door was opened for 30 s at 37°C	15 min	15 min
Electrical data		
Rated Voltage	230 V	120 V
Power frequency	50/60 Hz	60 Hz
Nominal power	0,3 kW	0,3 kW
Unit fuse	6,3 A	12,5 A
Phase (Nominal voltage)	1~	1~
Dimensions and weights		
Interior volume	62 L	62 L
Net weight of the unit (empty)	35 kg	35 kg
Load per rack	15 kg	15 kg
Permitted load	45 kg	45 kg
Wall clearance back	160 mm	160 mm
Wall clearance sidewise	100 mm	100 mm
Housing dimensions not incl. fittings and connections		
Width net	614 mm	614 mm
Height net	625 mm	625 mm
Depth net	577 mm	577 mm
Internal Dimensions		
Interior width	400 mm	400 mm
Interior height	440 mm	440 mm

Interior depth	350 mm	350 mm
Inner doors	1	1
Unit doors	1	1
<b>Environment-specific data</b>		
Energy consumption at 37°C	25 Wh/h	25 Wh/h
<b>Fixtures</b>		
Number of shelves (std./max.)	1/3	1/3

## DIMENSIONS Incl. fittings and connections [mm]

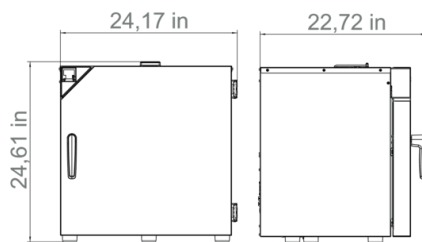

## ACCESSORIES

Designation	Description	BD-S 56	*	Art.-No.
Grid shelf	Incl. shelf supports			
	chrome plated	•	-	8012-1807
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	•	-	8012-2250
Rubber pads	set anti-slip feet	•	-	8012-2030

## SERVICES

Designation	Description	*	Art.-No.
<b>Calibration services</b>			
Temperature calibration	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate Travel to be quoted separate.	-	DL30-0110
	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate Travel to be quoted separate.	-	DL30-0102
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	-	DL30-0109
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate. Travel to be quoted separate. A total of 3 shelves are required.	-	DL30-0127

Designation	Description	*	Art.-No.
<b>Installation services</b>			
Unit installation	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work) Travel to be quoted as separate.	-	DL10-0100
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming. Travel to be quoted as separate Limited to 3 hours.	05	DL10-0500
<b>Maintenance services</b>			
BRONZE 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate. 10% Discount on Spare Parts and Labor for any additional Services needed Wear + tear parts to be purchased or billed to customer. Travel to be quoted as separate.	-	DL20-0720
GOLD 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate of Temperature. 25% Discount on Spare Parts and Labor for any additional Services needed. One additional courtesy visit per year can be used for service visit or operator. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	-	DL20-0915
Preventive maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Checking a test temperature specified by the user in center of usable space without certificate. Wear + tear parts need to be purchased or billed to customer. Travel to be quoted as separate.	-	DL20-0601
SILVER 3-year maintenance contract	Annual Device Inspection Includes Preventive Maintenance with Calibration Certificate 15% Discount on Spare Parts and Labor for any additional Services needed One additional courtesy visit per year can be used for service visit or operator training. Wear + tear parts to be purchased or billed to customer. Travel to be quoted separate.	-	DL20-0815
<b>Warranty service</b>			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	-	DL50-0010

**NOTES**

- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

**BINDER GmbH**

Tuttlingen, Germany  
TEL +49 7462 2005 0  
info@binder-world.com  
www.binder-world.com

**BINDER Inc.**

Bohemia, NY, USA  
TEL +1 631 224 4340  
usa@binder-world.com  
www.binder-world.us

**BINDER Environmental Testing  
Equipment (Shanghai) Co., Ltd.**

Shanghai, P.R. China  
TEL +86 21 685 808 25  
china@binder-world.com  
www.binder-world.com

**BINDER Asia Pacific (Hong Kong) Ltd.**

Kowloon, Hong Kong, P.R. China  
TEL +852 39070500  
asia@binder-world.com  
www.binder-world.com